


Tracker Alignment

Ready to perform the alignment procedure for the 3D tracker:

- Verify that the quasi-freefloating *Free-floating Restraint* is installed as shown with the *Chest Marker Support* attached with velcro to the restraint plate.

Note: Alignment of the marker structure on the restraint plate is critical. Please ensure that the top edge of the *Chest Marker Support* is aligned with the edge of the plate nearest the deck of the ISS module.

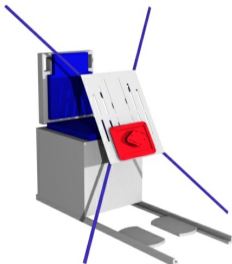

- Check that the field of view of the *Tracking Cameras* is free of obstructions.

Press '**Execute**' to perform alignment.

Action Completed Normally



Press '**Next**' to continue.

Action Completed with Errors



Press:

- **'Retry'** to repeat this action.
- **'Ignore'** to ignore and move on.
- **'Restart'** to restart this task or select another.

Recommendation: Attempt **'Retry'** no more than once or twice before consulting with POCC. **'Ignore'** and **'Restart'** should not be used without POCC authorization.