


## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the **GRIP Chair** is installed with the seat back deployed and the **Chest Marker Support**, **Hand Marker Support** and **Head Marker Support** are attached with velcro to the **GRIP Chair** as shown.

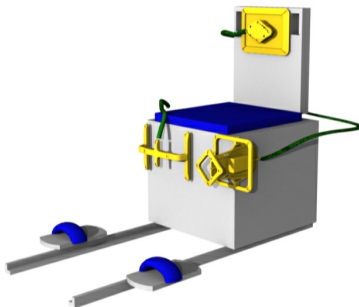

**Note:** Alignment of the **Chest Marker Support** on the backrest is critical. Please ensure that the top edge of the **Chest Marker Support** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the **Tracking Cameras**).

- Check that the field of view of the *Tracking Cameras* is free of obstructions.

Press '**Start**' to perform alignment.

# Tracker Alignment Running



**Follow the instructions on the screen.**

## HMD Assembly

- Attach the ***Head Marker Support*** to the front of the ***Virtual Reality Headset*** as shown:



- Secure with velcro straps on top and bottom.
- Attach the small junction box on the ***Head Marker Support*** cable to the top strap of the ***Virtual Reality Headset***.

**Note:** Alignment of the Marker Structure to the ***Virtual Reality Headset*** is critical. Please ensure that the structure is parallel to the front face of the ***Virtual Reality Headset***.

Press '**Next**' to continue.

## Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the ***GRIP Chair***.
- Temp stow the ***Virtual Reality Headset***, the ***Hand Marker Support*** and the ***Chest Marker Support*** to the right side of the ***GRIP Chair*** (subject's right side when seated facing the ***Tracking Cameras***).



- Ensure that the *Marker Array Cables* and *VR Headset Cable* pass behind the *GRIP Chair*.

Press '**Next**' to continue.