


## Tracker Alignment

Ready to perform the alignment procedure for the 3D tracker:

- Verify that the Quasi-Freefloating restraint is installed as shown with the Chest Marker Structure attached with velcro to the restraint plate.

**Note:** *Alignment of the marker structure on the restraining plate is critical. Please ensure that the top edge of the marker structure (see label) is aligned with the edge of the plate nearest the floor.*

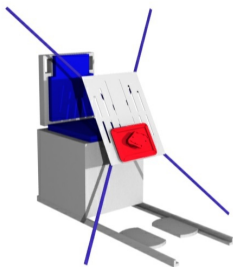

- Check that the field of view of the tracker cameras is free of obstructions.

Press '**Execute**' to perform alignment.

## Action Completed Normally



Press '**Next**' to continue.

## Action Completed with Errors



Press:

- **'Retry'** to repeat this action.
- **'Ignore'** to ignore and move on.
- **'Restart'** to restart this task or select another.

**Recommendation:** Attempt **'Retry'** no more than once or twice before consulting with POCC. **'Ignore'** and **'Restart'** should not be used without POCC authorization.