


## Tracker Alignment

Ready to perform the alignment procedure for the 3D tracker:

- Verify that the **GRIP Chair** is installed as shown with the seat back deployed and the **Chest Marker Support** attached with velcro to the backrest.

**Note:** Alignment of the **Chest Marker Support** on the backrest is critical. Please ensure that the top edge of the marker structure is aligned with the top edge of the backrest.

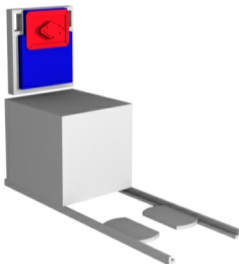

- Check that the field of view of the **Tracking Cameras** is free of obstructions.

Press '**Execute**' to perform alignment.

## Action Completed Normally



Press '**Next**' to continue.

## Action Completed with Errors



Press:

- **'Retry'** to repeat this action.
- **'Ignore'** to ignore and move on.
- **'Restart'** to restart this task or select another.

**Recommendation:** Attempt **'Retry'** no more than once or twice before consulting with POCC. **'Ignore'** and **'Restart'** should not be used without POCC authorization.