


## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the quasi-freefloating ***Floating Backrest*** is installed with the ***Chest Marker Support***, ***Hand Marker Support*** and ***Head Marker Support*** attached with velcro as shown to the restraint plate.

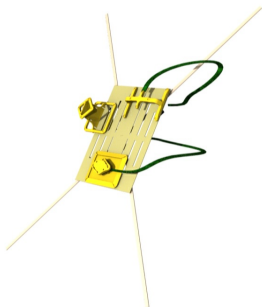

**Note:** Alignment of the ***Chest Marker Support*** on the restraint plate is critical. Check that the top edge of the ***Chest Marker Support*** is aligned with the edge of the plate nearest the deck of the ISS module and that the cable exits the support to your right when facing the plate.

- Check that the field of view of the ***Tracking Cameras*** is free of obstructions.

Press '**Start**' to perform alignment.

## Action Completed Normally



Press '**Next**' to continue.

## Action Completed with Errors



Press:

- **'Retry'** to repeat this action.
- **'Ignore'** to ignore and move on.
- **'Restart'** to restart this task or select another.

**Recommendation:** Attempt **'Retry'** no more than once or twice before consulting with COL-CC. **'Ignore'** and **'Restart'** should not be used without Col-CC authorization.