


## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the **GRIP Chair** is installed with the seat back deployed and the **Chest Marker Support**, **Hand Marker Support** and **Head Marker Support** are attached with velcro to the **GRIP Chair** as shown.

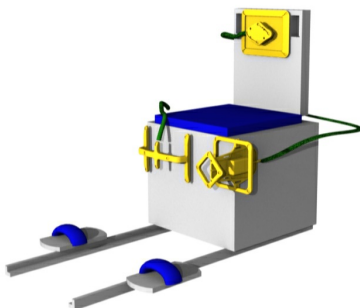


**Note:** Alignment of the **Chest Marker Support** on the backrest is critical. Please ensure that the top edge of the **Chest Marker Support** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the **Tracking Cameras**).

- Check that the field of view of the **Tracking Cameras** is free of obstructions.

Press '**Start**' to perform alignment.


# Tracker Alignment Running

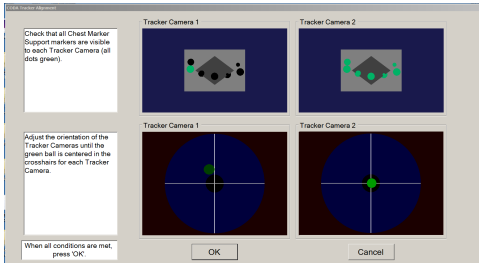

**Follow the instructions on the screen.**

## HMD Assembly

- Attach the ***Head Marker Support*** to the front of the ***Virtual Reality Headset*** as shown:

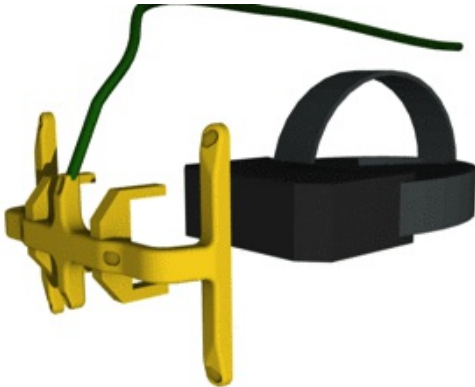

- Secure with velcro straps on top and bottom.
- Attach the small junction box on the ***Head Marker Support*** cable to the top strap of the ***Virtual Reality Headset***.


**Note:** Alignment of the Marker Structure to the ***Virtual Reality Headset*** is critical. Please ensure that the structure is parallel to the front face of the ***Virtual Reality Headset***.

Press '**Next**' to continue.

## Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the **GRIP Chair**.
- Temp stow the **Virtual Reality Headset**, the **Hand Marker Support** and the **Chest Marker Support** to the right side of the **GRIP Chair** (subject's right side when seated facing the **Tracking Cameras**).

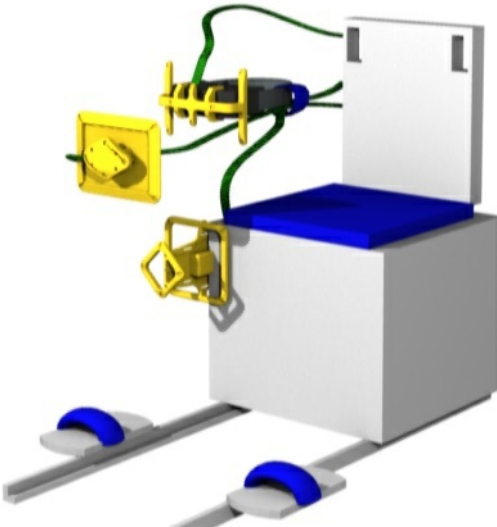

- Ensure that the **Marker Array Cables** and **VR Headset Cable** pass behind the **GRIP Chair**.

Press '**Next**' to continue.