


## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the quasi-freefloating ***Floating Backrest*** is installed with the ***Chest Marker Support***, ***Hand Marker Support*** and ***Head Marker Support*** attached with velcro as shown to the restraint plate.

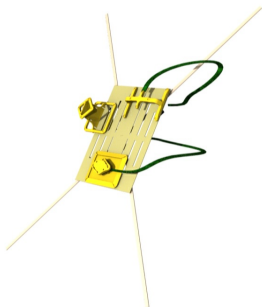


**Note:** Alignment of the ***Chest Marker Support*** on the restraint plate is critical. Check that the top edge of the ***Chest Marker Support*** is aligned with the edge of the plate nearest the deck of the ISS module and that the cable exits the support to your right when facing the plate.

- Check that the field of view of the ***Tracking Cameras*** is free of obstructions.

Press '**Start**' to perform alignment.


## Tracker Alignment Running

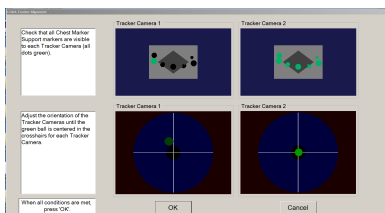

**Follow the instructions on the screen.**

## HMD Assembly

- Attach the *Head Marker Support* to the front of the *Virtual Reality Headset* as shown:

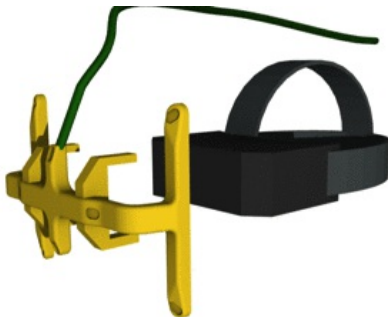

- Secure with velcro straps on top and bottom.
- Attach the small junction box on the *Head Marker Support* cable to the top strap of the *Virtual Reality Headset*.


**Note:** Alignment of the Marker Structure to the *Virtual Reality Headset* is critical. Please ensure that the structure is parallel to the front face of the *Virtual Reality Headset*.

Press 'Next' to continue.

## Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the **GRIP Chair**.
- Temp stow the **Virtual Reality Headset**, the **Hand Marker Support** and the **Chest Marker Support** to the right side of the **GRIP Chair** (subject's right side when seated facing the **Tracking Cameras**).

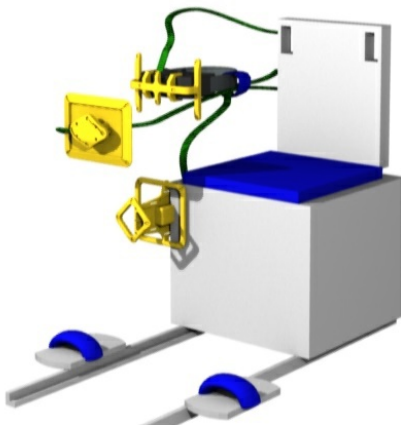

- Ensure that the **Marker Array Cables** and **VR Headset Cable** pass behind the **GRIP Chair**.

Press '**Next**' to continue.