


## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the **GRIP Chair** is installed with the seat back deployed and the **Chest Marker Support**, **Hand Marker Support** and **Head Marker Support** are attached with velcro to the **GRIP Chair** as shown.

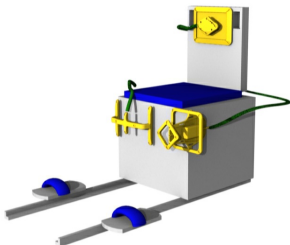

**Note:** Alignment of the **Chest Marker Support** on the backrest is critical. Please ensure that the top edge of the **Chest Marker Support** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the **Tracking Cameras**).

- Check that the field of view of the **Tracking Cameras** is free of obstructions.

Press '**Start**' to perform alignment.

## Action Completed Normally



Press '**Next**' to continue.

## Action Completed with Errors



Press:

- **'Retry'** to repeat this action.
- **'Ignore'** to ignore and move on.
- **'Restart'** to restart this task or select another.

**Recommendation:** Attempt **'Retry'** no more than once or twice before consulting with COL-CC. **'Ignore'** and **'Restart'** should not be used without Col-CC authorization.