


# GRASP Software User Interface


This document provides an overview of the steps followed by a subject or operator in order to run the GRASP installation and scientific protocols on the Perspectives laptop. It covers only the software aspects, assuming that the GRASP hardware (including shared GRIP hardware) is already installed per the PODF specific to those activities.

## *Starting GRASP on the Perspectives Laptop*

A crew member initiates Grasp activities by double clicking on the GRASP icon situated on the desktop. This will start up first a utility program to verify the installation of the virtual reality engine and VR Headset. The initial screen will be as follows while the VR system initializes:

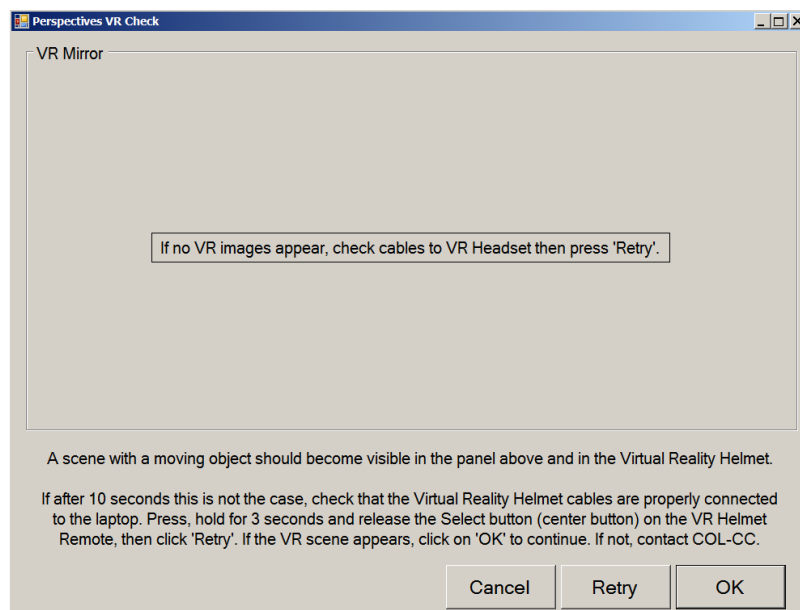


after 10-15 seconds a VR screen should appear as follows:

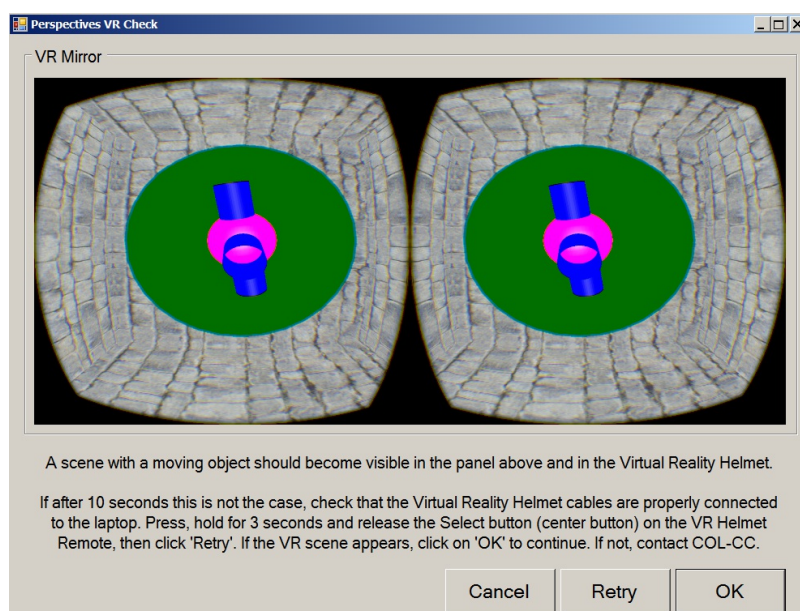


If the VR scene does not appear after 30 seconds, the operator should press and hold for 3 seconds the Select button on the VR Remote and then release, as indicated by the instructions on the screen. This should cause the VR scene to appear in the window and in the VR Headset. If not, the operator should contact Col-CC for further instructions.

After the VR engine has been started successfully, a window showing the status of the 3D tracker:

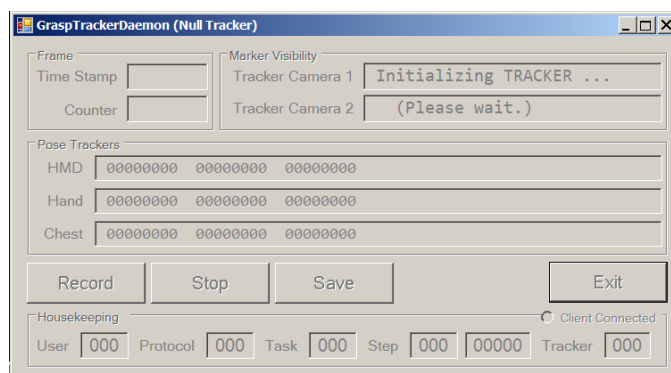


should appear for ~10 sec and then disappear (collapse to task bar).

Finally, the GRASP main menu should appear:

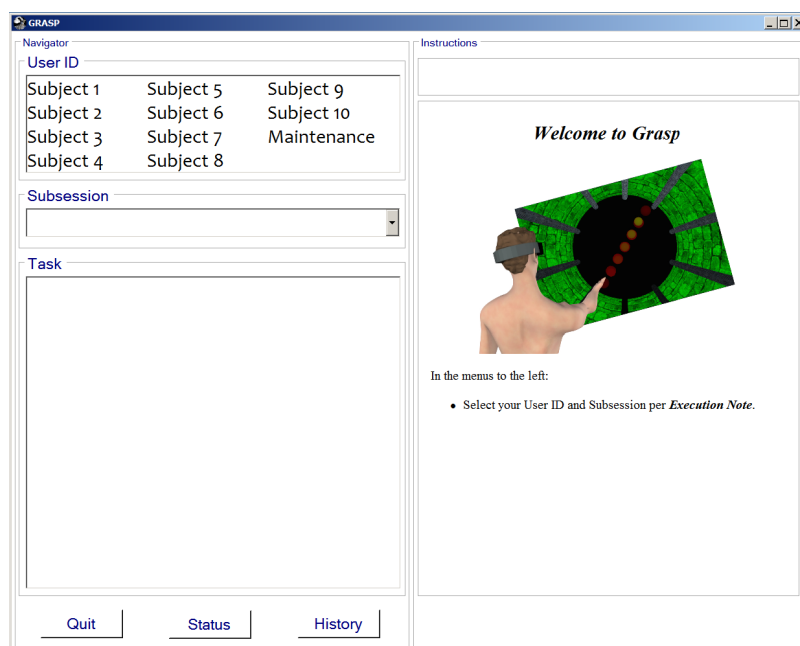


The subject or operator selects his or her designated subject ID in the upper-left panel of GRASP@ISS. Then the desired sub-session ID sub-session is selected via a pull-down menu:

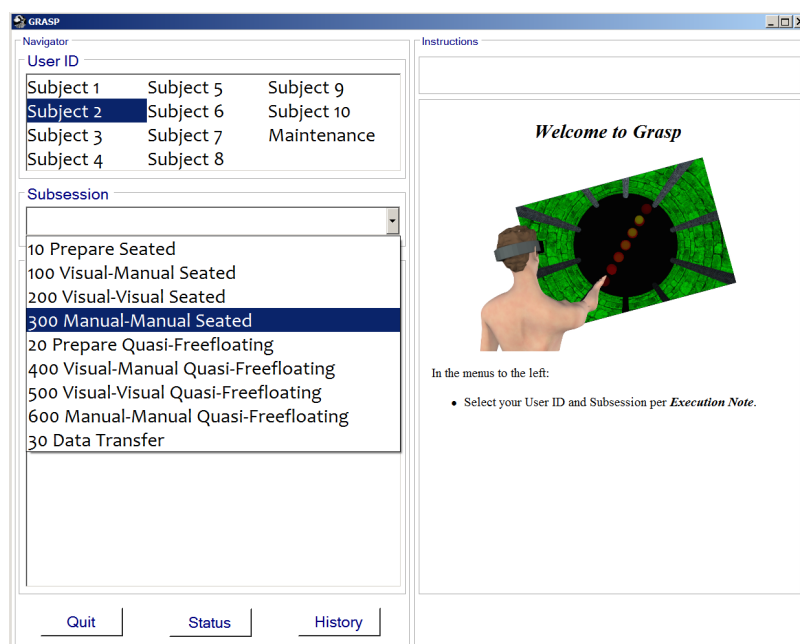


The subject ID and the sub-session to be performed at any given time is communicated to the crewmember via an Execution Note.

Selecting the sub-session protocol brings up a list of tasks to be executed in a defined order. The first task is automatically selected. An example is shown here:

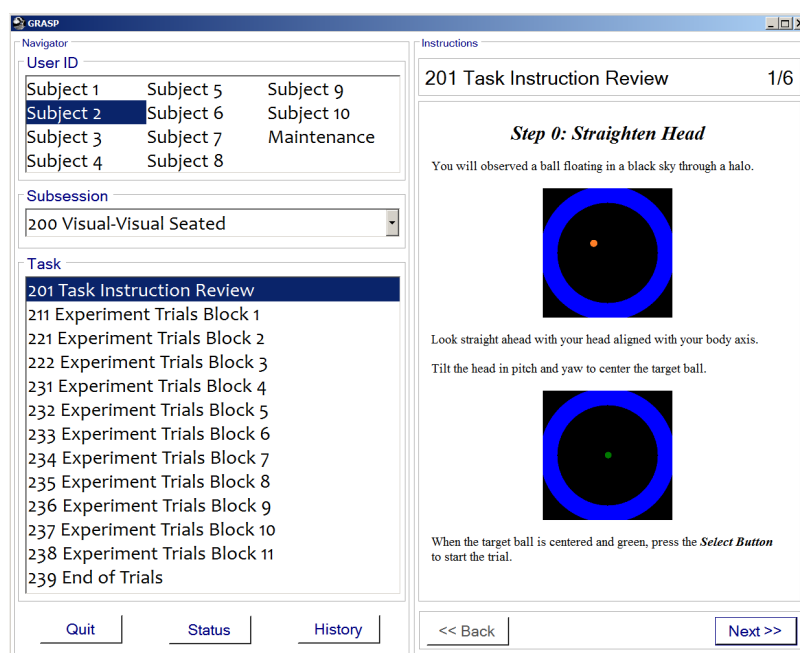

Note that at any time the subject can press the 'Status' button to bring up a screen showing the visibility of the markers on various marker structures and a test window showing the status of the VR hardware.

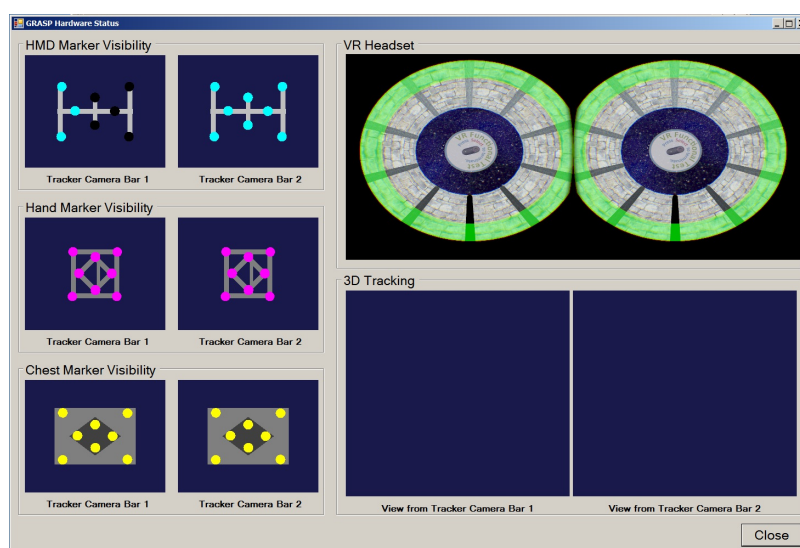


The subject can also press the 'History' button to display a list of the tasks that have been carried out on that day, with start and end time and an indication of the completion code for each task.

Start	Usr	Ssn	Tsk	Stp	Completion Status	End
2017.10.16 00:37:11	099	950	952	001	Completed normally (code 0)	2017.10.16 00:37:11
2017.10.16 00:39:06	099	950	952	001	Completed normally (code 0)	2017.10.16 00:39:06
2017.10.16 00:40:21	099	950	952	001	Completed normally (code 0)	2017.10.16 00:40:22
2017.10.16 00:41:03	099	950	952	001	Completed normally (code 0)	2017.10.16 00:41:03
2017.10.16 00:41:43	099	950	952	001	Completed normally (code 0)	2017.10.16 00:41:43
2017.10.16 00:46:03	099	950	951	001	Completed normally (code 0)	2017.10.16 00:46:03
2017.10.16 00:47:11	099	950	952	001	Completed normally (code 0)	2017.10.16 00:47:11
2017.10.16 01:31:42	099	950	954	001	Completed normally (code 0)	2017.10.16 01:31:42
2017.10.16 01:32:07	099	950	955	001	Completed normally (code 0)	2017.10.16 01:32:07
2017.10.16 01:34:16	099	950	954	001	Completed normally (code 0)	2017.10.16 01:34:16
2017.10.16 01:34:56	099	950	955	001	Completed normally (code 0)	2017.10.16 01:34:57
2017.10.16 01:35:44	099	950	952	001	Completed normally (code 0)	2017.10.16 01:35:44

## Executing Tasks

The panel on the right then guides the crew member through the steps required by each task. The panel may take one of four different forms.

1. First, the crew member is presented with one or more pages of instructions. Once the instructions have been read the subject proceeds to the subsequent page by pressing the 'Next' button. It is also possible to return to the previous step (except when at the first step) by pressing the 'Back' button.

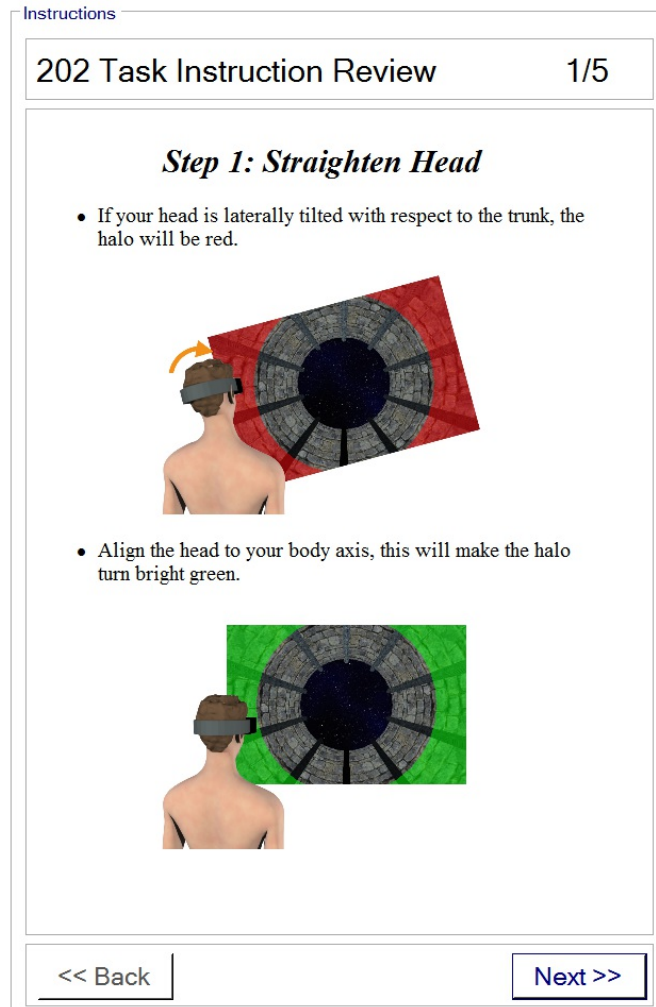

2. The crew member may then be prompted to launch a particular action. In the example shown here, the system is ready to perform the tracker alignment. The crew member presses the 'Execute' button to initiate the action. He or she may also press 'Back' to return to the previous action or instruction (except when it is the first step in a task). In exceptional cases, the subject may be told by COL-CC to press 'Skip' to move on to the next step without performing the action.

## Instructions

## 11 Preparation

1/4

**Tracker Alignment**

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the **GRIP Chair** is installed with the seat back deployed and the **Chest Marker Support**, **Hand Marker Support** and **Head Marker Support** are attached with velcro to the **GRIP Chair** as shown.



**Note:** Alignment of the **Chest Marker Support** on the backrest is critical. Please ensure that the top edge of the **Chest Marker Support** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the **Tracking Cameras**).

- Check that the field of view of the **Tracking Cameras** is free of obstructions.


Press '**Start**' to perform alignment.

&lt;&lt; Back

Skip &gt;&gt;

Start

3. Pressing 'Execute' typically causes a new program to execute, which will activate specific popup windows outside the GRASP@ISS main window. When execution of the task terminates, focus will return to the GRASP@ISS main window and the termination status of the task is displayed. If the task was successful, the crew member will typically see the following message allowing them to confirm and move on to the following step (by pressing 'Next') or to return and repeat the task (by pressing 'Back').

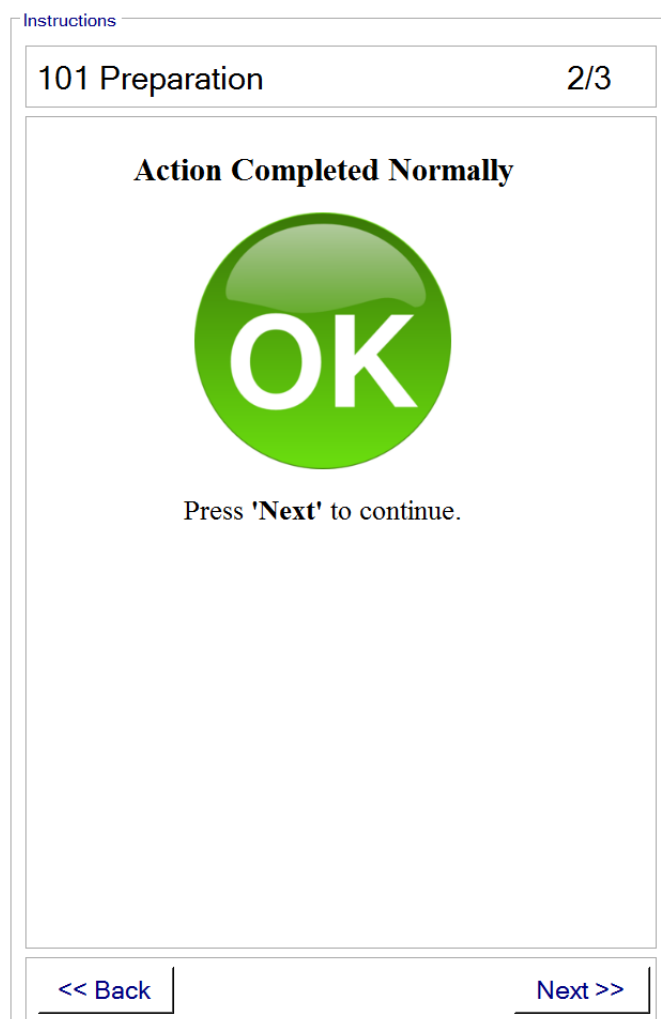

4. If the task did not complete successfully, an alert will be displayed that shows the error return code (in the 'Error Code' box) and that typically includes a graphical or text explanation of the anomaly. Shown here is the generic message that may be used for any anomalous condition, but more specific messages may also be displayed (see below).

Instructions

203 Practice Trials1/1

**Action Completed with Errors**  


Press:

- **'Retry'** to repeat this action.
- **'Ignore'** to ignore and move on.
- **'Restart'** to restart this task or select another.

**Recommendation:** Attempt **'Retry'** no more than once or twice before consulting with POCC. **'Ignore'** and **'Restart'** should not be used without POCC authorization.

Error Code  [Restart](#) [Ignore](#) [Retry](#)

In these circumstances, the crew member is invited to repeat the task by pressing 'Retry'. If a second attempt also fails, the crew member has the option to move on to the next step anyway, by pressing 'Ignore' or to repeat the entire task from the beginning by pressing 'Restart'.

The full set of possible instruction, status and alert screens are included in the following sections.