


Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the quasi-freefloating *Floating Backrest* is installed with the *Chest Marker Support*, *Hand Marker Support* and *Head Marker Support* attached with velcro as shown to the restraint plate.

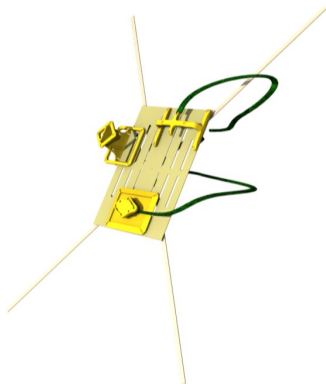


Note: Alignment of the *Chest Marker Support* on the restraint plate is critical. Check that the top edge of the *Chest Marker Support* is aligned with the edge of the plate nearest the deck of the ISS module and that the cable exits the support to your right when facing the plate.

- Check that the field of view of the *Tracking Cameras* is free of obstructions.

Press '**Start**' to perform alignment.


Tracker Alignment Running

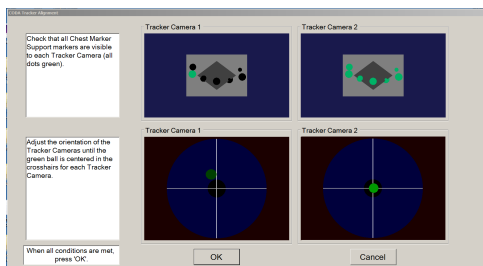

Follow the instructions on the screen.

HMD Assembly

- Attach the ***Head Marker Support*** to the front of the ***Virtual Reality Headset*** as shown:

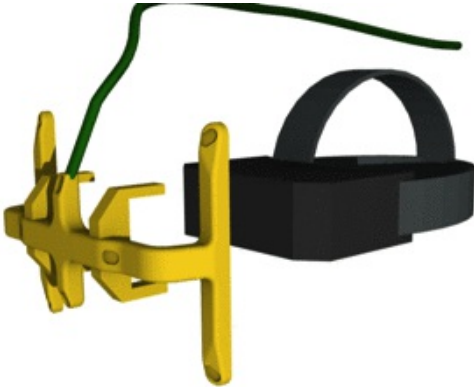

- Secure with velcro straps on top and bottom.
- Attach the small junction box on the ***Head Marker Support*** cable to the top strap of the ***Virtual Reality Headset***.


Note: Alignment of the Marker Structure to the ***Virtual Reality Headset*** is critical. Please ensure that the structure is parallel to the front face of the ***Virtual Reality Headset***.

Press '**Next**' to continue.

Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the **GRIP Chair**.
- Temp stow the **Virtual Reality Headset**, the **Hand Marker Support** and the **Chest Marker Support** to the right side of the **GRIP Chair** (subject's right side when seated facing the **Tracking Cameras**).

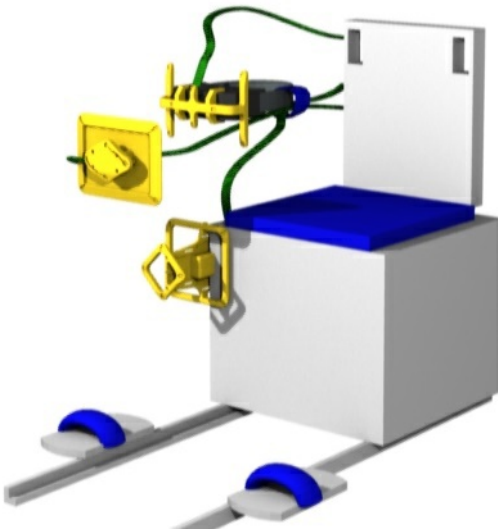

- Ensure that the **Marker Array Cables** and **VR Headset Cable** pass behind the **GRIP Chair**.

Press '**Next**' to continue.